




# PRESSALIT®

<b>Product</b> <b>Item number</b>	PLUS shower seat 450 with aperture, manually height and manually sideways adjustable R7455112
<b>Colour</b>	VVS: 777516412 NRF: 6317136 EAN: 5708590359984 HMI: 124242 FDA: INN/D135865 Anthracite grey
<b>Product description</b>	 <p>PLUS shower seats are offered in a number of colours. They are available as both fixed height, electrically adjustable and manually adjustable, with and without aperture, and with an option of back and armrests. They are suitable for wall-mounting or for fitting to Pressalit PLUS wall tracks.</p> <p>PLUS shower seat 450 with aperture, with backrest. Manually height adjustable by 380 mm and sideways adjustable, without the use of tools. The seat is foldable. Materials: Comfortable polyurethane foam over a plywood core, fitted to an aluminium profile. For fitting to horizontal PLUS wall track (not included).  <i>Replaces R7441</i></p>
<b>Adjustability</b> <b>Adjustment height</b> <b>Product capacity</b>	Height and sideways adjustable 380 mm Maximum user mass: 200 kg. When the seat is loaded with 200 kg, each wall fastener is loaded with 23 kg. <i>The seat has been tested statically with safety factor 1.5, which is a requirement according to the applied assistive products standards.</i>
<b>Materials</b> <b>Ergonomy</b>	<p>The shower seat is made of polyurethane (PU) foam with a plywood core and this is covered with an in-mould coated surface that is protected using a clear PU lacquer. The backrest is made of polyurethane (PU) foam protected using a clear PU lacquer. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS.</p> <p>Ergonomically designed seat - contact surface are 450 x 437 mm. Soft and slip-resistant material that has a warm feel.</p>

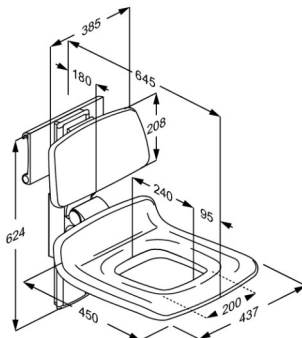

# PRESSALIT®

<b>Product</b>	PLUS shower seat 450 with aperture, manually height and manually sideways adjustable
<b>Item number</b>	R7455112
<b>Features</b>	The seat can fold up when not in use or to maximise floor space. The seat and backrest can be adjusted simultaneously thanks to a safe and easy click system.
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
<b>Cleaning</b>	<p>PLUS shower seat.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p><b>DECALCIFICATION OF ALUMINUM PARTS:</b> Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p>
<b>Declaration of conformity</b>	CE MDR This product is in compliance with all relevant requirements of: - Regulation (EU) 2017/745 on medical devices.

[pressalit.com](http://pressalit.com)

# PRESSALIT®