




# PRESSALIT®

<b>Product</b> <b>Item number</b>	PLUS backrest, height and sideways adjustable R6444035
<b>Colour</b>	NRF:6317051 EAN:5708590365527 HMI: 124205  Red (RAL 3003)
<b>Product description</b>	 <p>PLUS backrests are offered in a number of colours and versions, they are available as both fixed height and adjustable. They are suitable for wall-mounting behind a toilet without cistern or for mounting on one of Pressalit's toilet lifters.</p> <p>PLUS backrest, 385 x 209 mm, soft feel.          Height adjustable 150 mm and sideways adjustable, for mounting on a horizontal walltrack behind a wall-mounted toilet without cistern.          Can be removed easily when not needed and reattached again without the use of tools.          Materials: Aluminium and polyurethane foam with a plywood core.          For fitting to horizontal PLUS wall track (not included).</p>
<b>Adjustability</b> <b>Adjustment height</b>  <b>Product capacity</b>	Height and sideways adjustable  150 mm  Maximum load: 50 kg. Maximum user weight: 100 kg. When the backrest is loaded with 50 kg, each wall fastener will be carrying 33 kg. <i>Backrests need to be able to support the weight of the person who is leaning against it.</i> <b><i>The backrest is not intended to carry the full body weight of a person, and the test methods applicable to backrests reflect this: ISO 17966 prescribes a static test at 50% of maximum user weight, i.e. when the maximum load on a wash basin bracket is given as 50 kg, it is – according to ISO 17966 – the equivalent of a user weighing 100 kg leaning on the backrest.</i></b>
<b>Materials</b>  <b>Ergonomy</b>	The backrest is made of polyurethane (PU) foam which is covered with an in-mould coated surface that is protected using a clear PU lacquer.  Soft and slip-resistant material that has a warm feel.

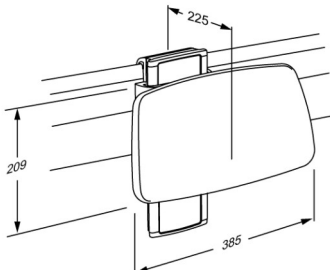

# PRESSALIT®

<b>Product</b>	PLUS backrest, height and sideways adjustable
<b>Item number</b>	R6444035
<b>Features</b>	For toilet without cistern
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.
<b>Cleaning</b>	<p>PLUS backrest.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p><b>DISINFECTANTS:</b></p> <p>Erisan OXY+</p> <p>Tvättspit 70%</p> <p>CaviCide surface disinfectant decontaminant cleaner</p> <p>Rodalon Desinfektionsspray</p> <p>Terralin protect (SCHÜLKE+)</p> <p>Dismonzon Plus (BODE Chemie GmbH)</p> <p><b>NON-DISINFECTANTS:</b></p> <p>Wetrok Caledor (Wetrok)</p> <p>Petrosol Citro (Büfa)</p> <p>Important: Polyurethane surfaces of PLUS shower seats and PLUS backrests contains an antibacterial surface, which will be compromised if cleaned with anything other than soap and water.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p><b>DECALCIFICATION OF ALUMINUM PARTS:</b></p> <p>Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p> <p>CE MDR This product is in compliance with all relevant requirements of: - Regulation (EU) 2017/745 on medical devices.</p> <p><b>Expected service life:</b> 10 years.</p> <p><b>Contraindications for use:</b> It is not recommended that people who need to lie down while being changed use the backrest for support.</p>
<b>Declaration of conformity</b>	



# PRESSALIT®

<b>Product</b>	PLUS backrest, height and sideways adjustable
<b>Item number</b>	R6444035

---

[pressalit.com](http://pressalit.com)

---

Pressalit A/S  
Pressalitvej 1  
DK-8680 Ry  
Tel.: +45 8788 8788  
Fax: +45 8788 8789  
E-mail: [pressalit@pressalit.com](mailto:pressalit@pressalit.com)

# PRESSALIT®