




PRESSALIT®

Product Item number	PLUS fold down grab bar with integrated counter-balance, 27.6", right hand operated R371470108
Color	NRF:6191115 EAN:5708590368429 HMI: 123566 FDA: ILS/D135866 Dark blue (RAL 5011)
Product description	 <p>All PLUS fold down grab bars come with an integrated counter-balance, which ensures a controlled movement between horizontal and vertical positions. Suitable for use with toilet pans projecting 21.3/21.7". <i>Replaces discontinued models R3106, RF106, R3127 og RF127.</i></p> <p><i>PLUS fold down grab bar, counter-balanced. Projection 24.8", ergonomically shaped to provide a clear sense of where the fold down grab bar ends. Right hand operated. Height adjustable without using tools. Detachable and foldable.</i></p> <p><i>Materials: The arms are constructed of 10-micrometer anodized aluminium, with a non-slip gripping surface made of soft SEBS incorporating a glass-reinforced nylon core</i></p> <p><i>For mounting on vertical track or SELECT TL2 side profile (not included).</i></p>
Adjustability Adjustment height Product capacity	Height adjustable 9,46 inch Maximum load: 330 lbs. When the fold down grab bar is loaded with 330 lbs, each wall fastener in the vertical wall track is loaded with 90 lbs. <i>PLUS fold down grab bars have been tested statically with 1540 lbs (close to the wall), 935 lbs (in the middle) and 495 lbs (the outermost point).</i>
Materials	The grab bars are primarily constructed in 10-micrometer anodized aluminium, with a non-slip gripping surface at the end made of soft SEBS and incorporating a glass-reinforced nylon core. The hinge arm is also made of soft SEBS. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS polymer.

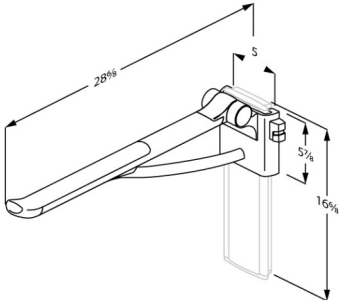

PRESSALIT®

Product	PLUS fold down grab bar with integrated counter-balance, 27.6", right hand operated
Item number	R371470108
Ergonomy	The oval shape gives an easy and comfortable surface to lean and push against. The end of the support arm is enlarged so that it is easier to pull forward from a sitting position and also provides a clear sense of where the support arm ends. Soft and slip-resistant material that has a warm feel.
Features	Folds against the wall. The underside is shaped for the mounting of various accessories.
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date.



PRESSALIT®

Product	PLUS fold down grab bar with integrated counter-balance, 27.6", right hand operated
Item number	R371470108
Cleaning	<p>This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Oxyplus-acid 2% (Diversey) WipeClean UniversalCleaning Chlorine Disinfection Wipes</p> <p>NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS: Pressalit recommends that water deposits (Calcium, lime etc.) be removed from wall tracks and other aluminum parts using a solution of three parts of 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts should be wiped over with a damp cloth wrung out in plain water.</p>
Declaration of conformity	This product is in compliance with all relevant requirements of regulation (EU) 2017/745 on medical devices.



PRESSALIT®

Product	PLUS fold down grab bar with integrated counter-balance, 27.6", right hand operated
Item number	R371470108
	<p>Intended use The fold down grab bar is used as a support for sitting down or getting up from a toilet, or for sideways transfer from a wheelchair.</p> <p>Operating environment For indoor use in a wet environment, such as a bathroom or a washroom, in private homes, nursing homes and hospitals. For institutional use, the expected life time is 10 years. <i>Life cycle test performed according to ISO 17966:2016.</i></p> <p>For maintenance, please refer to operation and maintenance manual, which is enclosed with the product, and which can also be found on the website under "Downloads".</p>

pressalit.com

Pressalit A/S
Pressalitvej 1
DK-8680 Ry
Tel.: +45 8788 8788
Fax: +45 8788 8789
E-mail: pressalit@pressalit.com

PRESSALIT®