




PRESSALIT®

| | |
|--------------------------------------|--|
| Product Item number | PLUS support arm with integrated counter-balance, 850 mm, with operating panel, left hand operated R372585112 |
| Colour | Anthracite grey |
| Product description | <p>VVS:615878212 NRF:6191143 EAN:5708590368689 HMI: 138174</p>  <p>All PLUS support arms come with an integrated counter-balance, which ensures a controlled movement between horizontal and vertical positions.</p> <p>PLUS support arm with operating panel, counter-balanced, 850 mm. Height adjustable 150 mm. Foldable. Left handed. Only for use with TL3/TL4 toilet lifter equipped with toilet pan projecting 540/550 mm. Ergonomically shaped to provide a clear sense of where the support arm ends. Materials: The arms are constructed of 10-micrometre anodized aluminium, with a non-slip gripping surface made of soft SEBS incorporating a glass-reinforced nylon core. Excl. mounting kit.</p> |
| Adjustability | Height adjustable |
| Adjustment height | 150 mm |
| Product capacity | Maximum load: 150 kg |
| Materials | <p>The arms are primarily constructed in 10-micrometre anodized aluminium, with a non-slip gripping surface at the end made of soft SEBS and incorporating a glass-reinforced nylon core. The hinge arm is also made of soft SEBS. Internal components and fasteners are primarily made of stainless steel. The hinge barrel is capped with ABS polymer.</p> |
| Ergonomy | <p>The oval shape gives an easy and comfortable surface to lean and push against.</p> <p>The end of the support arm is enlarged so that it is easier to pull forward from a sitting position and also provides a clear sense of where the support arm ends.</p> <p>Soft and slip-resistant material that has a warm feel.</p> |

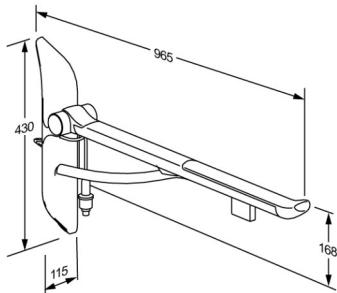

PRESSALIT®

| | |
|--------------------|--|
| Product | PLUS support arm with integrated counter-balance, 850 mm, with operating panel, left hand operated |
| Item number | R372585112 |
| Features | The underside is shaped for the mounting of various accessories. |
| Warranty | Pressalit makes a commitment to remedy material and manufacturing defects for a period of five years from the invoicing date. |
| Cleaning | <p>PLUS support arms.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container.</p> <p>DISINFECTANTS:</p> <p>Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Oxyplus-acid 2% (Diversey) WipeClean UniversalCleaning Chlorine Disinfection Wipes</p> <p>NON-DISINFECTANTS:</p> <p>Suma Star free D1 Taski 100 Taski 200 Milizid Citro</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p> <p>DECALCIFICATION OF ALUMINUM PARTS:</p> <p>Pressalit recommends that chalky deposits be removed from wall tracks and other aluminum parts using a solution of three parts 30% acetic acid to seven parts tap water. When the aluminum parts have been washed down with this solution, the parts must be wiped over with a damp cloth wrung out in clear water.</p> |



PRESSALIT®

| | |
|---|---|
| Product Item number | PLUS support arm with integrated counter-balance, 850 mm, with operating panel, left hand operated R372585112 |
| WEEE Declaration of conformity | <p>Please enquire of your local Pressalit subsidiary or distributor (in countries where there is no subsidiary) that will arrange for recycling of electronic waste in accordance with the WEEE directive.</p> <p>CE MDR-RoHS This product is in compliance with all relevant requirements of:</p> <ul style="list-style-type: none">- Regulation (EU) 2017/745 on medical devices- Council Directive 2011/65/EC on RoHS. <p>Intended use The support arm is used as a support for sitting down or getting up from a toilet, or for sideways transfer from a wheelchair.</p> <p>Operating environment For indoor use in a wet environment, such as a bathroom or a washroom, in private homes, nursing homes and hospitals. For institutional use, the expected life time is 10 years. <i>Life cycle test performed according to ISO 17966:2016.</i></p> <p>For maintenance, please refer to user guide, which is enclosed with the product, and which can also be found on the website under "Downloads".</p> |

pressalit.com