




# PRESSALIT®

<b>Product</b> <b>Item number</b>	MCT 2 change trolley, electrically height adjustable R9432721000
<b>Colour</b>	VVS: 777542220 EAN: 5708590373218 HMI: 128432  White
<b>Product description</b>	 <p>Pressalit offers a variety of mobile and wall-hung, height-adjustable and fixed height shower and changing tables of high quality. They are designed to provide a stable platform for showering and changing adults and children with reduced mobility. These tables can be used in schools, institutions, hospitals and private homes and are all designed to be durable, comfortable and easy to clean. The design and flexibility of the tables ensures optimal working conditions for the carer as well as a high level of comfort and safety for the user.</p> <p>MCT 2 change trolley, 820 x 2190 mm, two fold-down sides.  Electrically height adjustable (battery operated) with wired remote control.  400 mm lifting range, from 550 to 950 mm.  Equipped with central lock system, wire basket for convenient storage, bumpers on all four corners as well as battery and charger.  Lifting columns and wired remote control IPX6, battery and charger IPX5.  Mattress to be ordered separately.</p>
<b>Adjustability</b> <b>Adjustment height</b>  <b>Product capacity</b>  <b>Speed upwards</b>	Height adjustable  400 mm  Maximum user weight (MUM): 200 kg. Maximum load (SWL): 250 kg. (user weight + water + accessories). <i>This product has been tested statically with a load of 625 kg as the IEC 60601-1 standard requires testing with safety factor 2.5.</i>  Speed up mm: 19 mm/s Speed down mm: 19 mm/s
<b>Materials</b>	Aluminium, plastics, compact laminate and powder-lacquered stainless steel.

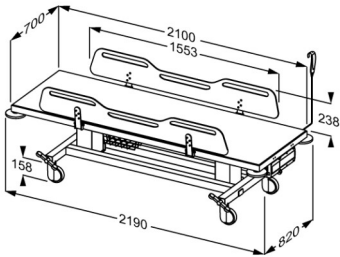

<b>Product</b>	MCT 2 change trolley, electrically height adjustable
<b>Item number</b>	R9432721000
<b>Ergonomy</b>	<ul style="list-style-type: none"><li>• The ergonomically designed safety rails have a practical lower central section, which provides better access and work environment for the carer while transferring, changing and dressing the user.</li><li>• One hand operated safety rails allows the carer to maintain eye contact as well as physical contact with the user during the whole process, the foldable safety rails also provide easy access from all four sides of the table.</li><li>• The cut-out grips in the safety rails enable users to assist when being turned, which motivates physical training.</li><li>• A storage drawer enables the carer to store supplies and place what is needed for the shower immediately at hand in advance.</li><li>• Easy to clean, smooth contact surfaces which do not collect dirt and bacteria.</li></ul>
<b>Features</b>	<ul style="list-style-type: none"><li>• The trolley can be moved from room to room for easy transport and multiple users. Two lifting columns are used to create an aesthetically pleasing, hygienically optimized, physically and visually stable construction.</li><li>• Intuitive, easy-to-use drop down rails with one hand operation.</li><li>• The centralized brake system allows the carer to simultaneously lock all four wheels by stepping on any single activator, ensuring better safety and ease of use in any position.</li><li>• All locking and release activators are activated by stepping down on them from above to ensure ergonomic interaction, and to prevent foot and toe injuries. Non-slip surface and pattern on all brakes prevent slipping even in wet conditions.</li><li>• A rotating bumper on each of the trolley's corners protects surrounding walls and the product itself when in contact and during relocation.</li><li>• Powered height adjustment with battery and charger. The wall-hung charger for the battery can be mounted in any room, incl. outside the bathroom if so desired. When the battery is running low, a beep will sound each time the wired remote control is activated.</li><li>• Battery charging time 3-4 hours.</li><li>• Turning diameter: 2310 mm.</li></ul>
<b>Mains connection (V)</b>	230
<b>Actuator (motor)</b>	BL1 (Linak) IPX6
<b>Control box</b>	COBO / CU20 (Linak) IPX6
<b>Wired remote control</b>	HB30 (Linak) IPX6
<b>Power consumption (maximum load)</b>	Power consumption by maximum load, watt: 100 W



# PRESSALIT®

<b>Product</b>	MCT 2 change trolley, electrically height adjustable
<b>Item number</b>	R9432721000
<b>Warranty</b>	Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is placed in an atmosphere heavy in salt or chloride, the warranty will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product.
<b>Cleaning</b>	<p>General.</p> <p>The product is to be cleaned with hot water or a mild soapy cleaning agent. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged.</p> <p>The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.</p>
<b>Declaration of conformity</b>	<p>CE MDR-MASK-RoHS This product is in compliance with all relevant requirements of:</p> <ul style="list-style-type: none"><li>- Regulation (EU) 2017/745 on medical devices</li><li>- Council Directive 2006/42/EC on machinery</li><li>- Council Directive 2011/65/EC on RoHS.</li></ul>

[pressalit.com](http://pressalit.com)

# PRESSALIT®