




PRESSALIT®

Product Item number	SCT 1000 shower change table, electrically height adjustable R8593118000
Color	EAN:5708590377667 FDA:INQ/D124319 White
Product description	 Pressalit offers a variety of mobile and wall-hung, height-adjustable and fixed height shower and changing tables of high quality. They are designed to provide a stable platform for showering and changing adults and children with reduced mobility. These tables can be used in schools, institutions, hospitals and private homes and are all designed to be durable, comfortable and easy to clean. The design and flexibility of the tables ensures optimal working conditions for the carer as well as a high level of comfort and safety for the user. SCT 1000 shower change table, powered, height adjustable 15.8"; from 18.9" - 34.7", with wired hand control Foldable. Incl. safety rail. Outer measurements 29.9" x 70.9". Materials: Rilsan coated steel and PVC-coated canvas (phthalate free). Excl. mounting kit.
Adjustability Adjustment height Product capacity	Height adjustable 15,76 inch Maximum user weight (MUM): 441 lbs. Maximum load (SWL): 441 lbs. When the table is loaded with 441 lbs, each mounting bolt supports 60 lbs. <i>This product has been tested statically with a load of 830 lbs as the IEC 60601-1 standard requires testing with a safety factor 2.5.</i>
Speed up Noise level	Speed up: 0,43 inch/s Speed down: 0,35 inch/s 50 dB
Materials	The table is made from practical and durable PVC canvas which is phthalate free. The frame is powder-coated steel. Material: Steel, PVC.

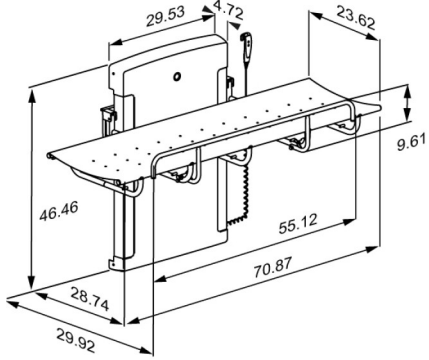

Product	SCT 1000 shower change table, electrically height adjustable
Item number	R8593118000
Ergonomy	Caregivers can easily set the height of the table, enabling them to obtain the correct working position. This ensures safe showering and changing for the user as well as minimizing the risk of back injury for the caregiver.
Features	Electric with wired hand control – the red buttons clearly indicate “Up” and “Down”. The wired hand control comes with a holder for fitting to the wall. The integrated pneumatic cylinders enable the table to be folded out of the way when not in use to maximize floor space in front of the unit. Splash curtain (optional) for the changing table protects the caregiver against splashing.
Mains connection (V)	120
Actuator (motor)	LA40 Careline (Linak), 4000 N, IPX6
Control box	CA63 Careline (Linak), IPX6
Wired remote control	HB30 (Linak), IP66
Power consumption (maximum load)	Power consumption by maximum load: 390 W Power consumption by standby: 0,5 W
Warranty	Pressalit makes a commitment to remedy material and manufacturing defects for a period of three years from the invoicing date. NOTICE If the product is placed in an atmosphere heavy in salt or chloride, the warranty will be reduced from the normal 3 years to 1 year. Salt and chloride have a corrosive effect on the products components and will cause a decrease in the life expectancy of the product.



PRESSALIT®

Product	SCT 1000 shower change table, electrically height adjustable
Item number	R8593118000
Cleaning	This product should be cleaned with hot water or a mild cleaning solution. Do not use chlorine based, corrosive or abrasive cleaners as it may result in the surface being damaged. The following cleaning agents have been tested by Pressalit and may be used, as long as they are diluted according to the instructions on the container. DISINFECTANTS: Clinimax Difficil-S Erisan OXY+ Sprit (husholdningssprit) Suma Bac D10 Virkon Rodalon Desinfektionsspray Lysoformin Spezial CaviCide Sufanios Premium Tvättsprit 70% Accel PREvention Wipes Lingettes NON-DISINFECTANTS: Suma Star free D1 Taski 100 Taski 200 Milizid Citro The product must be wiped down with a soft rag wrung out in clear water immediately after cleaning, as some materials may otherwise be discolored and stained, and the surface coating may deteriorate.
Declaration of conformity	This product is solely intended for the US and Canada, and has no EU declaration of conformity.



Product	SCT 1000 shower change table, electrically height adjustable
Item number	R8593118000

Mounting in different wall types

	Concrete >B25 V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.	
	Solid Brick >= Mz12 V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.	
	Solid Sand Lime Brick >= Hlz12 V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. To be ordered separately.	
	Aerated Concrete P4 V8658 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs GB14. To be ordered separately.	
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with wood V8608 - Mounting kit, 6 lag bolts Ø10 x 70 - plugs SX12. NB! Plugs are not needed. To be ordered separately.	



PRESSALIT®

Product Item number	SCT 1000 shower change table, electrically height adjustable R8593118000		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with 22 mm plywood Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		
	Drywall (2 x 12,5 mm) water-resistant, re-enforced with steel plates (.08 inch) Wall fasteners (6 pcs) must be bought separately, for drywall it is possible to use e.g. SnapTight or WingIt anchors or similar.		
	Drywall (2 x 12,5 mm) water-resistant Mounting in this wall type is not recommended - risk of fatal injury to persons and walls.		

It is the responsibility of the installer to ensure that suitable wall fasteners are used.

pressalit.com